


Saline Liquid Cal


Lee Drake
June 18, 2015




Magnesium

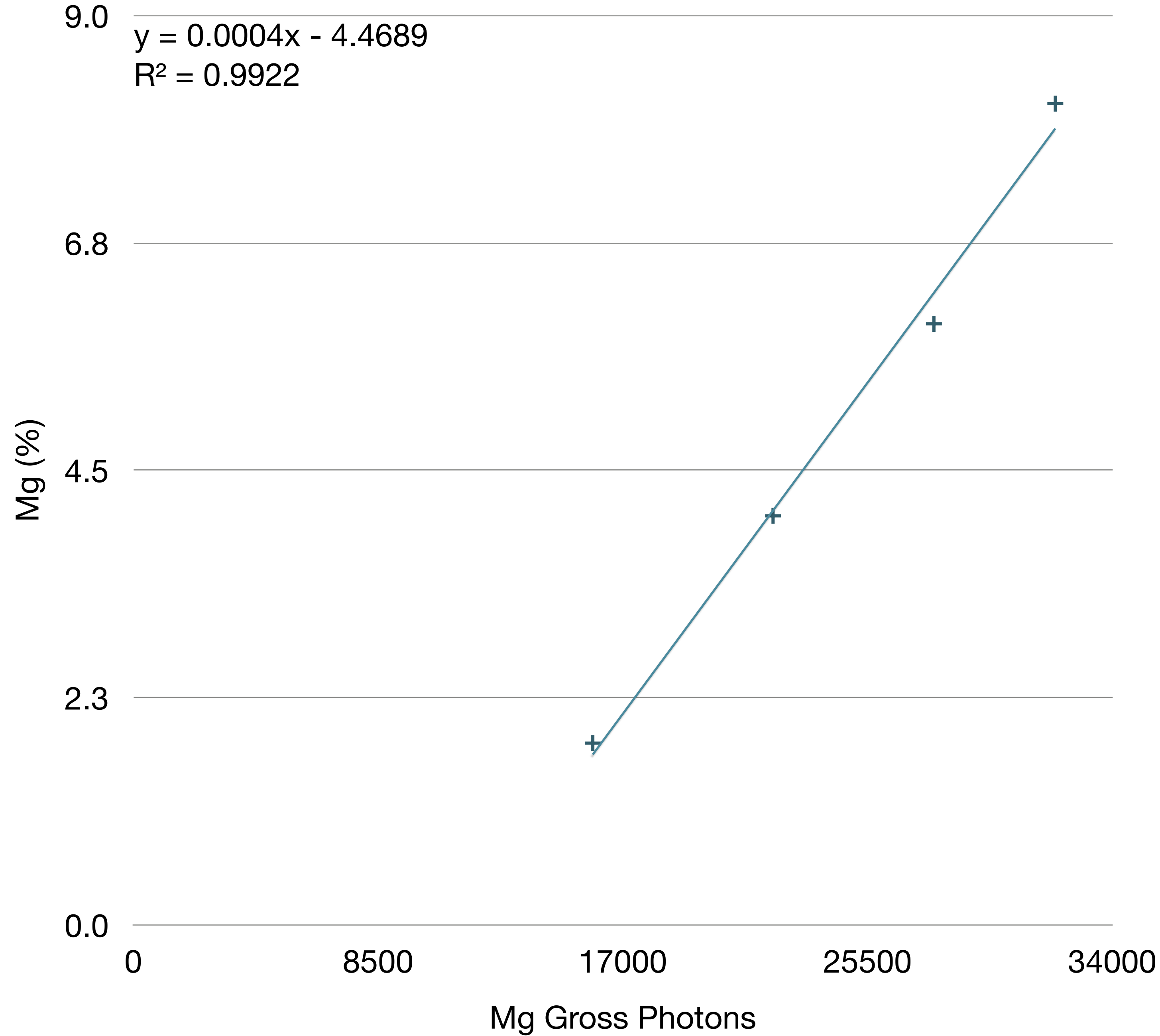


Sulfur

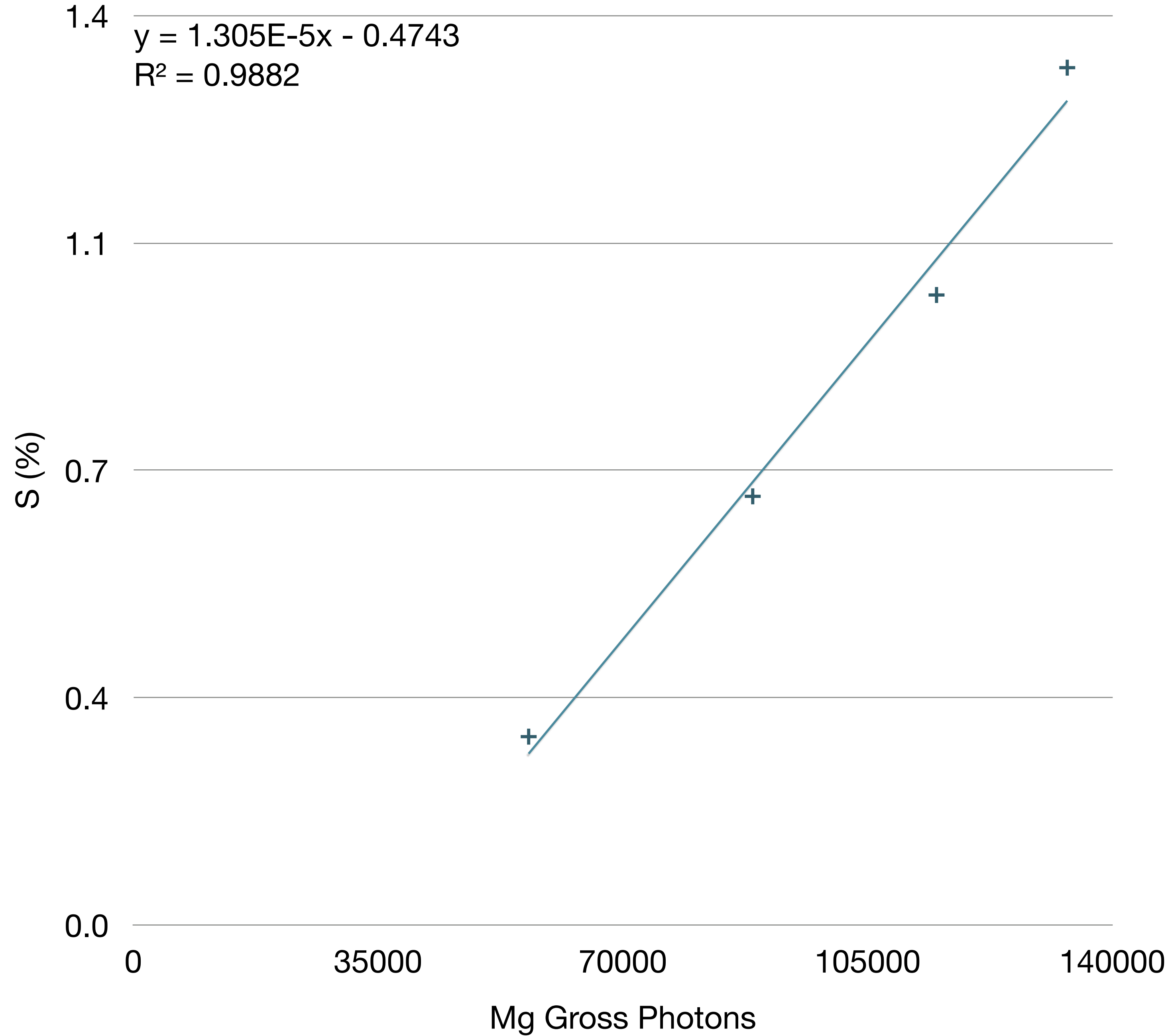


Chlorine

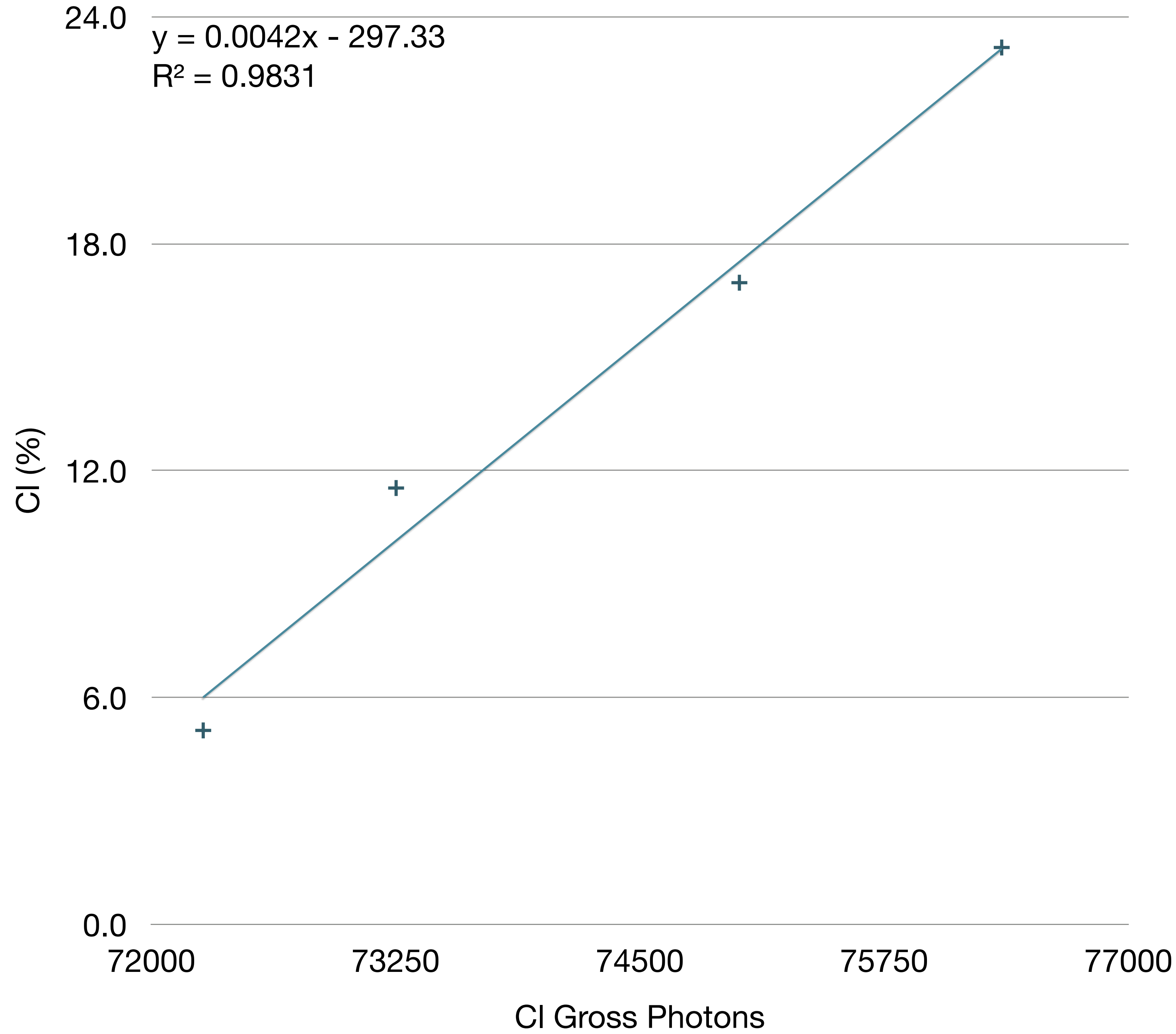

Potassium

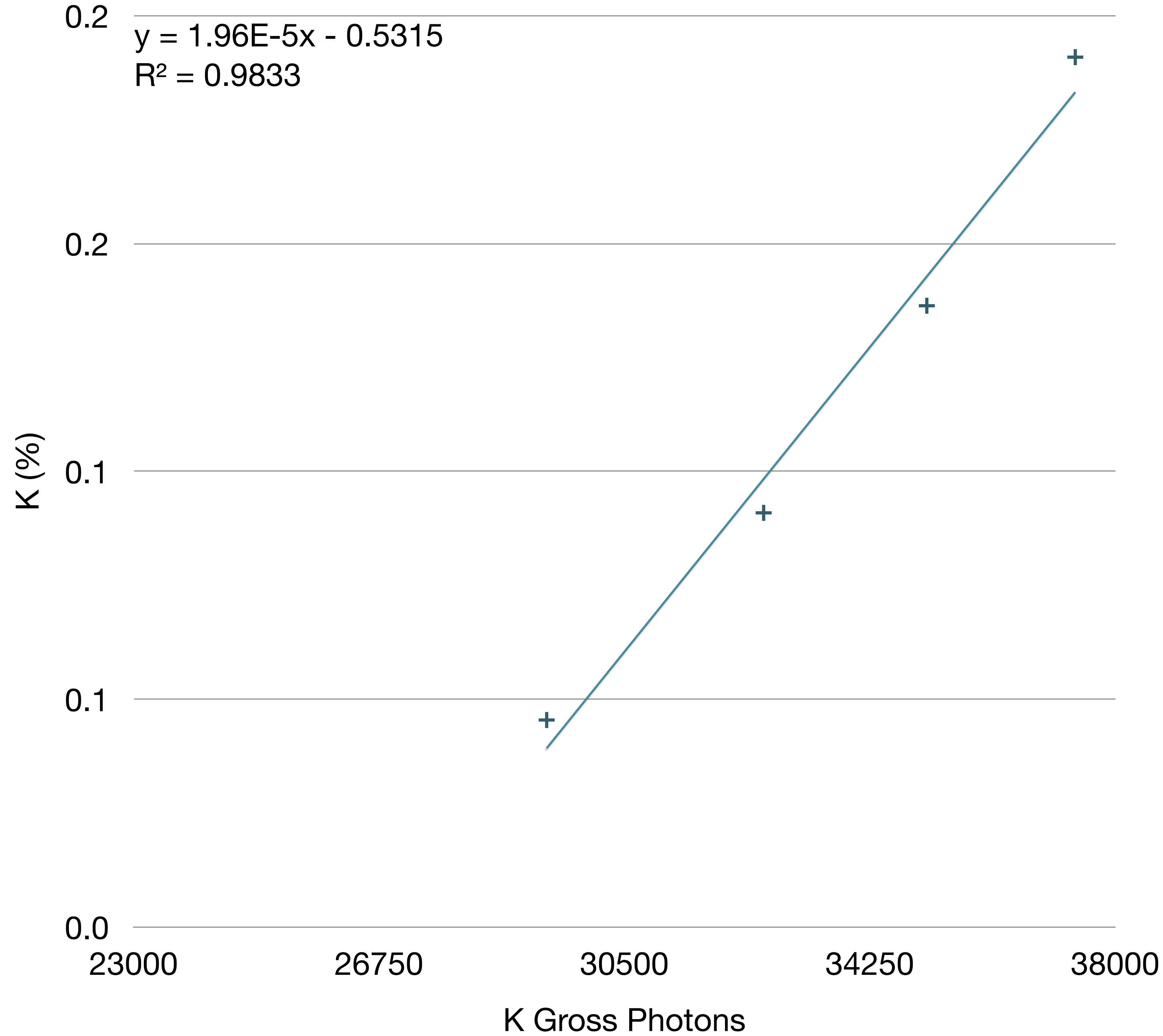

Sodium

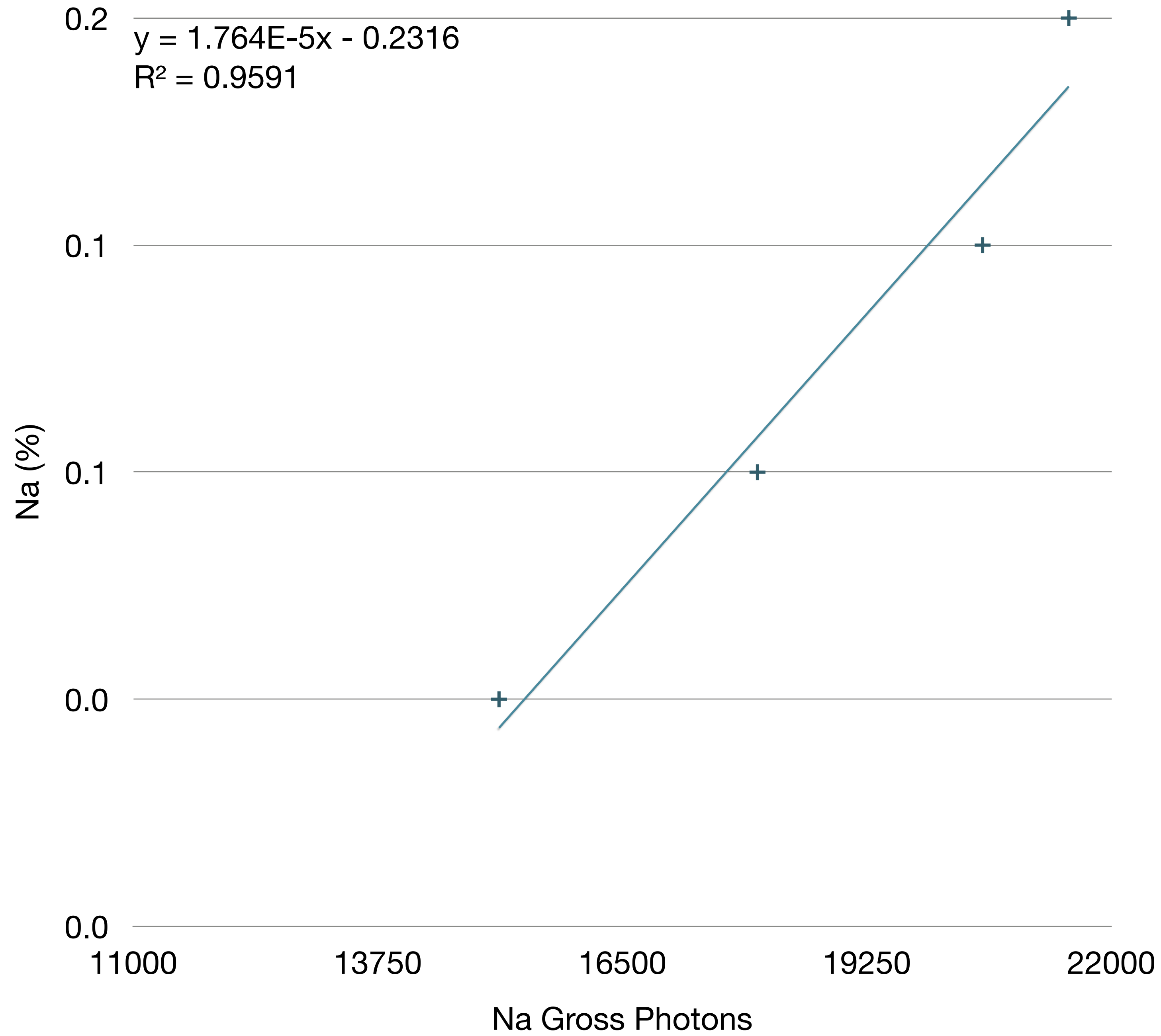

Given Values

Na	Mg	S	Cl	K
0.16	31.84	1.32	23.19	0.21
0.12	23.31	0.97	16.97	0.15
0.08	15.86	0.66	11.54	0.10
0.04	7.05	0.29	5.13	0.05

Quantified XRF Results

NaKa1	MgKa1	S Ka1	ClKa1	K Ka1
0.16	32.44	1.34	23.46	0.21
0.12	22.77	0.95	16.76	0.15
0.08	14.90	0.62	11.09	0.10
0.04	7.95	0.32	5.52	0.05

